

SIEMENS

syngo MammoReport

SP

System Manual

Wiring Diagrams

© Siemens AG 2005

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Register 6

Print No.: SPB7-250.844.30.01.02

English

Doc. Gen. Date: 10.05

Chapter	Page	Revision
All	All	01

Document revision level

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

Copyright notice

The DICOM functionalities in the *syngo* MammoReport system are based on the libraries from Merge Technologies Inc.

Celsius is a trademark of Fujitsu Siemens Computers.

Windows XP Professional is a trademark of Microsoft Corporation.

syngo is a trademark of Siemens AG.

Coronis and MediCal Pro are trademarks of Barco.

Dome and CXtra are trademarks of Planar Systems, Inc.

Softcopy reading software by MeVis BreastCare,
copyright © 2001 - 2005 by MeVis BreastCare.



	Page
1 System overview	1-1
General	1-1
syngo MammoReport system overview	1-1
2 Cables	2-1
Cable connections.	2-1
3 Connectors	3-1
Workstation connectors (CRT Monitors / for R610 only)	3-1
Workstation connectors (Siemens/Barco/Planar TFT Monitors)	3-2
4 Changes to previous version	4-1


Page

General

This document describes the main functions of the *syngo* MammoReport system and the electrical connections between the parts.

syngo MammoReport system overview

An overview of the *syngo* MammoReport system is provided below.

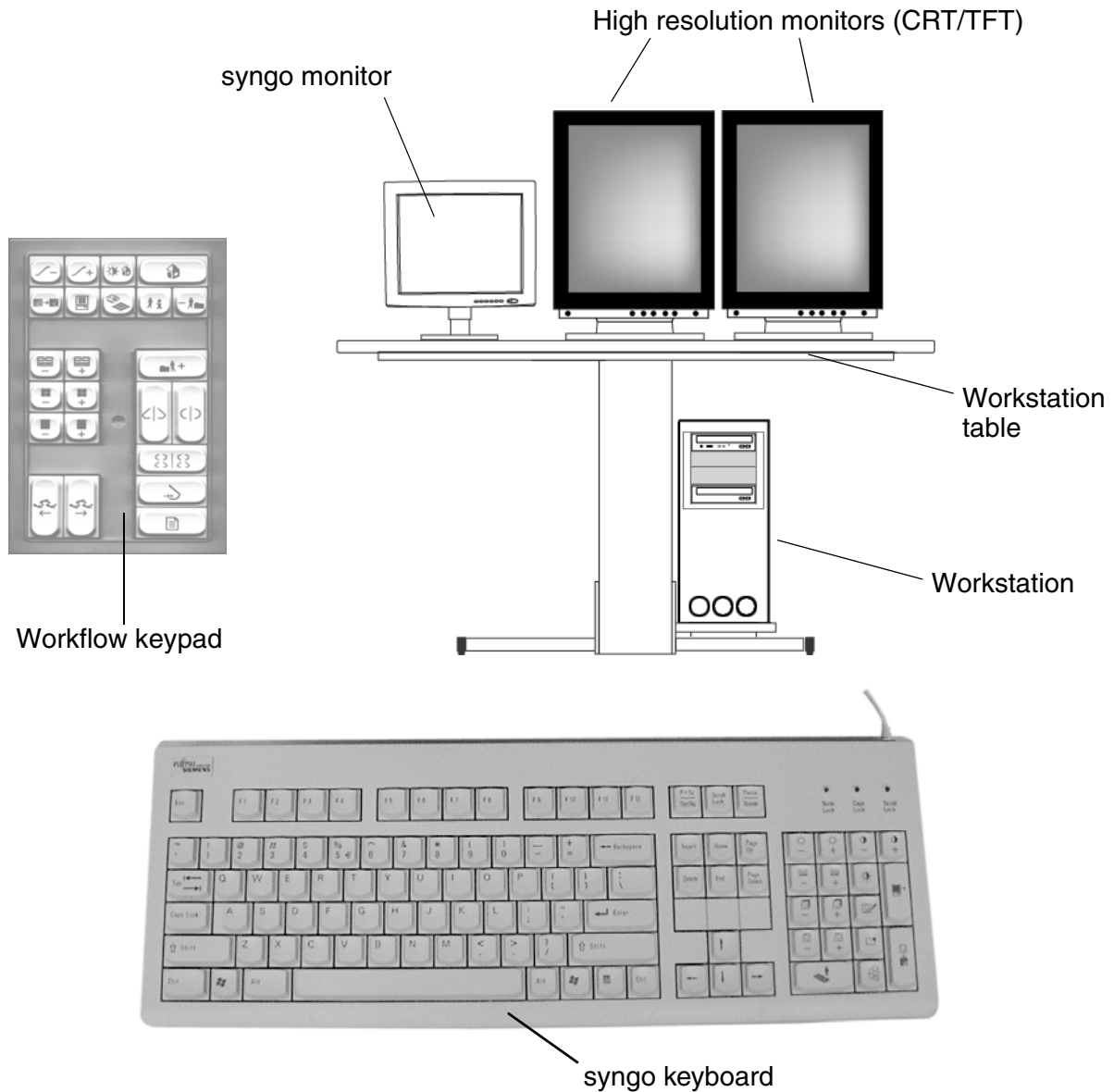

The *syngo* MammoReport system's major parts are:

- 2 Displays (either CRT or TFT)
 - CRTdisplays: SIMOMED SMM21201P
 - TFT displays: DSB 2103-D-5MP, Barco MFGD 5621 HD , or Planar (model C5i)
- Celsius workstation
- Workflow keypad
- *syngo* MammoReport software
- syngo keyboard
- syngo monitor
- Windows XP professional compatible printer, provided by the customer

Additional options to the *syngo* MammoReport system are:

- One additional 147 GB hard disk (R610)
- One additional 300 GB hard disk (R630)

NOTE

For actual system options, see price book.

Cable connections

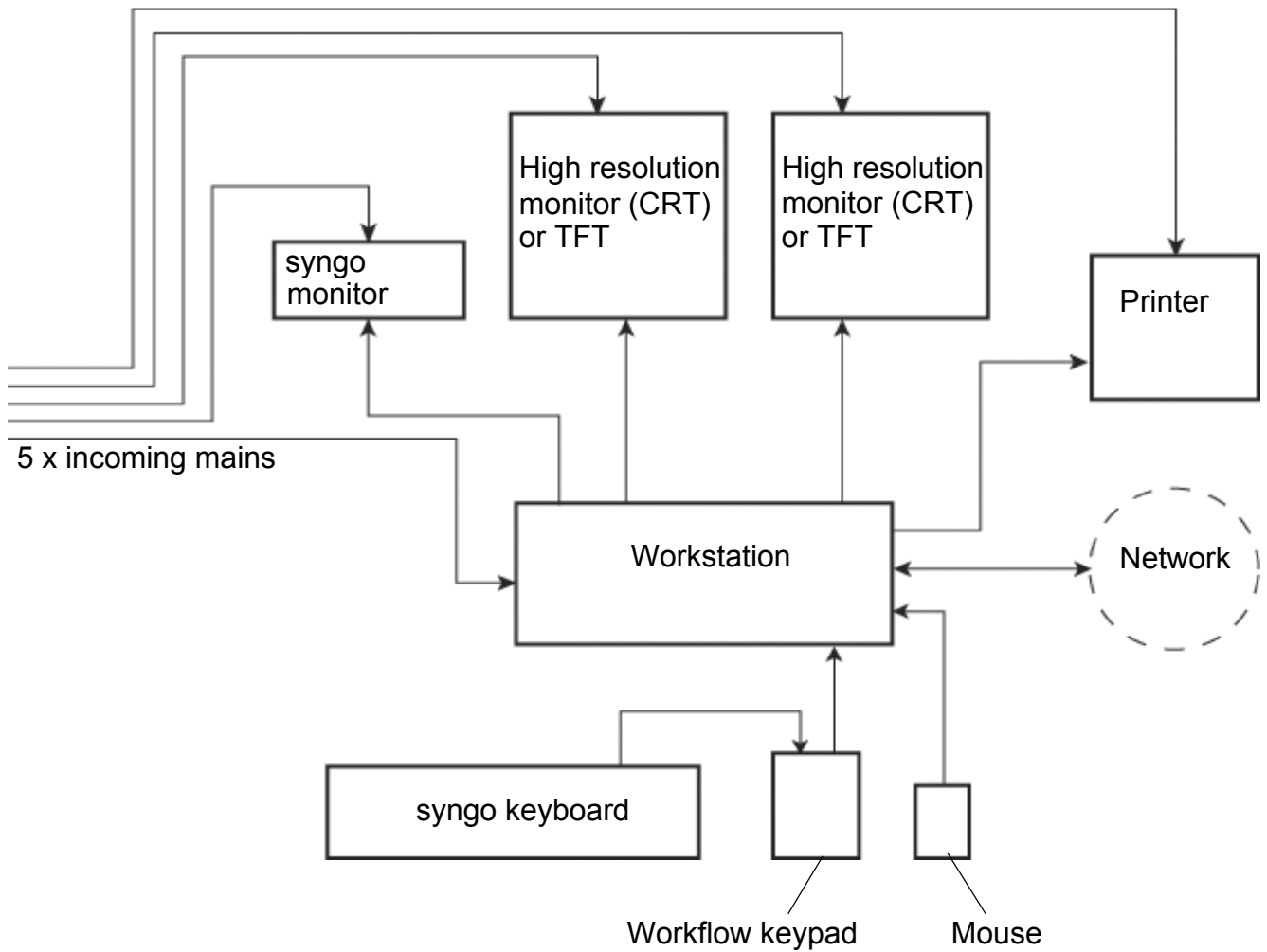


Fig. 1 Cable connections

NOTE

The *syngo* Mammography system can be connected to line voltage 100-240 V.

This page intentionally left blank.

Workstation connectors (CRT Monitors / for R610 only)

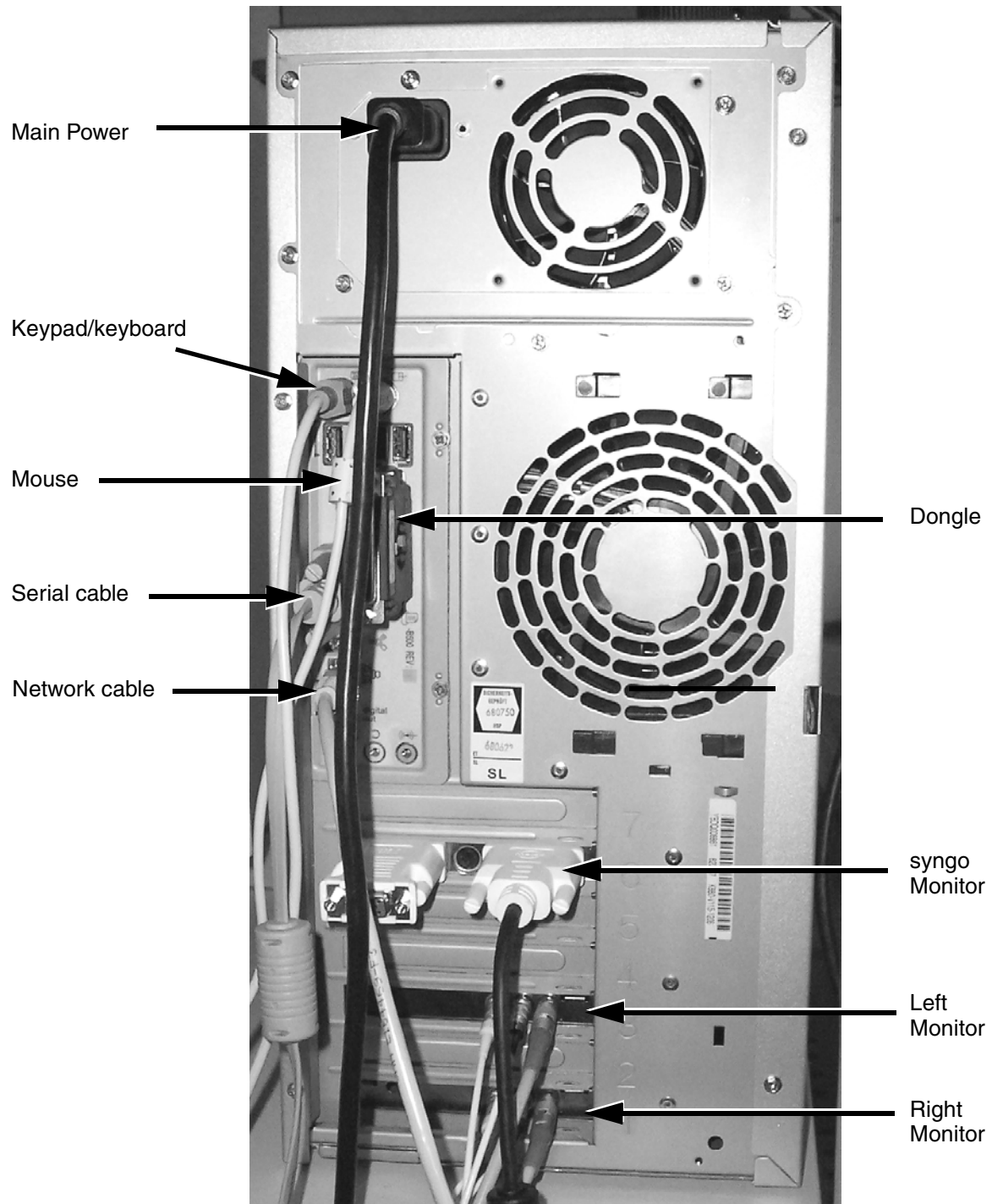


Fig. 1 Connectors Celsius (R610)

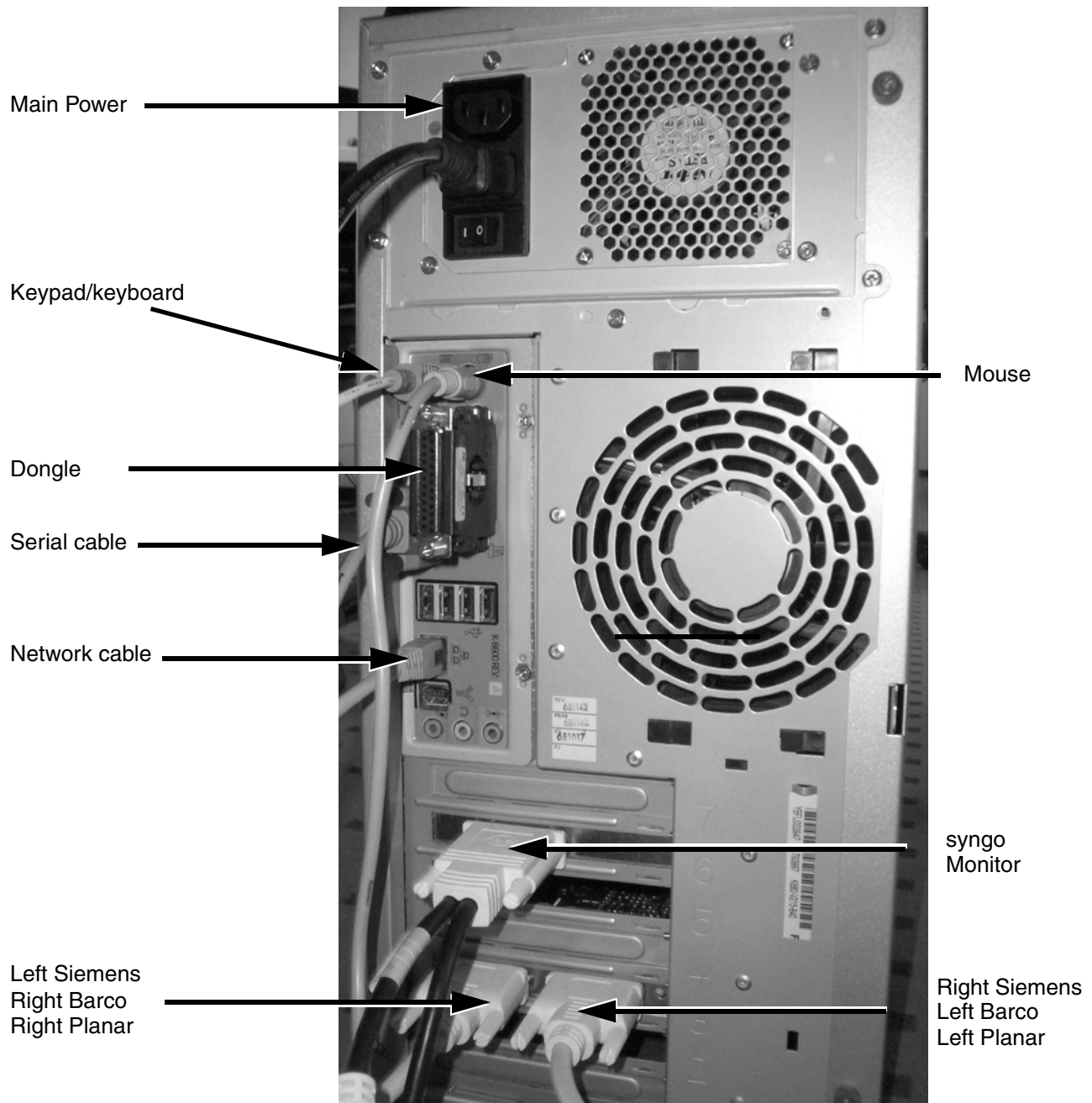
Workstation connectors (Siemens/Barco/Planar TFT Monitors)

Fig. 2 Connectors Celsius R630

Chapter	Page	Change
all		initial version

This page intentionally left blank.